 TABLE DE LIT TASMANIE

 TASMANIE BEDSIDE TABLE

 MESITA DE CAMA TASMANIA

REF. 823021

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INSTRUCTIONS MANUAL EN

IMPORTANT

Read the instructions carefully before using this product and keep them if necessary.

User: If you are unable to read or understand any warnings, precautions or instructions, please contact your health care professional or distributor before using this product or risk and damage may result.


Consult your doctor or health care professional before using a walking aid. He/she will advise and inform you on the most appropriate walking technique for you.

Distributor: this leaflet must be given to the users of this product.

PRODUCT DESTINATION

Tasmanie bedside tables are medical devices. They allow the user to read or take a meal in bed or in an armchair. They are adjustable in height and tilt to adapt to the needs of the user during the day.

PRODUCT DESCRIPTION

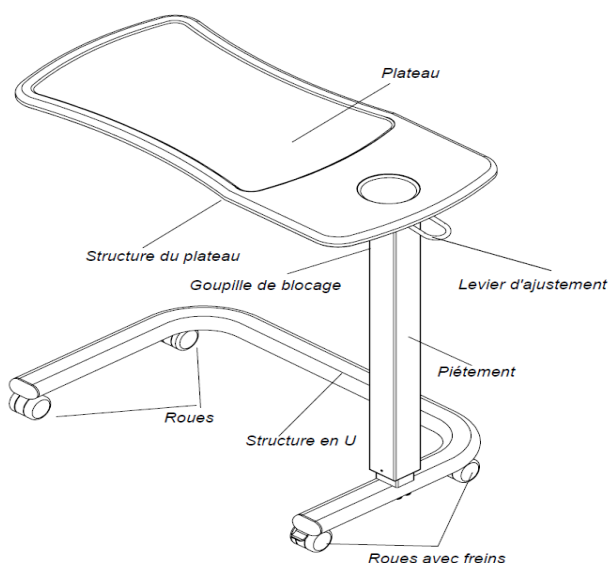Galerie de plateau : **Tray Gallery**

Grand plateau basculant : **Large tilting tray**

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SAFETY PRECAUTIONS


Failure to observe the precautions for use may result in the risk of serious injury.

- Use the bedside table only as described in this manual.
- Install the bed table on a flat surface.
- Store the bedside table in an indoor environment, away from moisture and dust.
- This bedside table uses a high-pressure pneumatic cylinder to adjust the height of the bed top. Do not remove the locking pin from the column assembly before use.
- Before use, make sure that all parts are present and in good condition. In case of damage, please contact your dealer.
- Carefully follow the assembly instructions and ensure that all components are properly secured before use.
- If any parts are damaged, broken or worn, contact your dealer.

MOUNTING

INSTALLATION

Place the bedside table on a clear surface and remove the packaging. Do not discard the packaging before the assembly is completed.

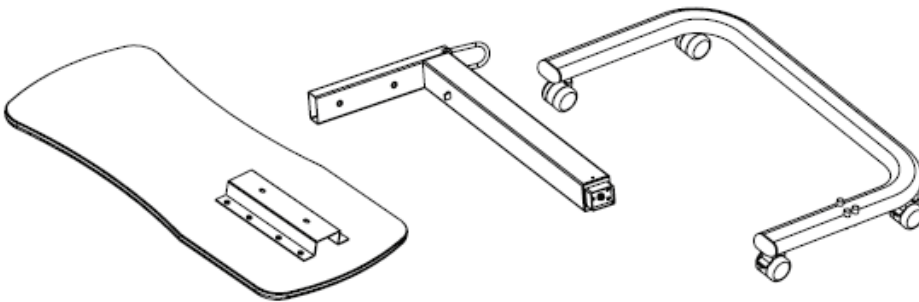
STEP 1


Place the table upside down on a flat surface as shown.

Insert the base above the table top frame, aligning the screw holes as shown

Insert the hexagon socket screws through the frame into the table top frame as indicated.

Use the M8 hexagonal key to tighten the screws.


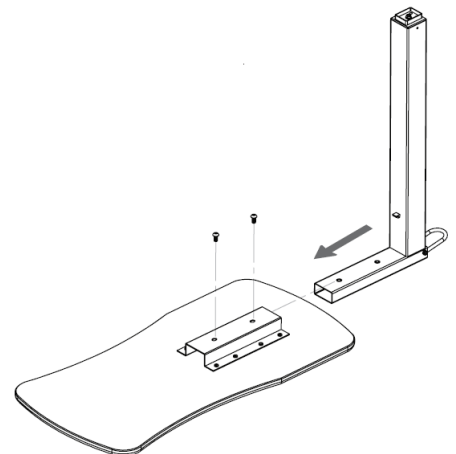
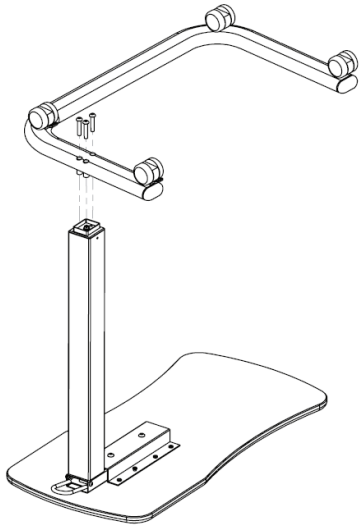
STEP 2

Place the U-shaped frame upside down over the base as indicated on the left. Screw the 3 screws through the plate, the base and the end of the base.

STEP 3

Turn the table over and place it on the floor. Make sure that all parts are perfectly tightened before using the bed table. Remove the locking pin from the actuator.

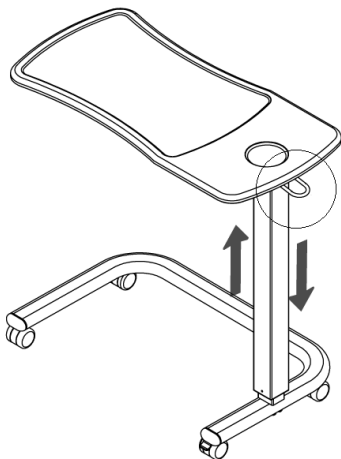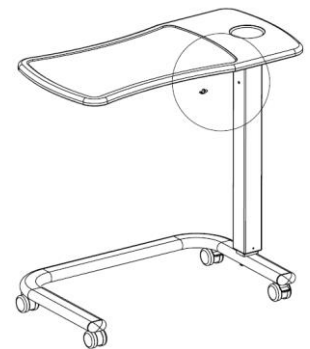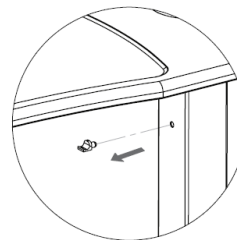HEIGHT ADJUSTMENT

Simply pull the table top upwards to adjust the table height. To lower the table, pull the adjustment lever, adjust to the correct height and then release the lever.

CARE AND MAINTENANCE

Failure to care for and maintain your product may result in the risk of serious injury.

- Clean the product regularly with warm water and a soft cloth. Do not use abrasive products, detergents or solvent-based products.
- Do not use cleaning products containing phenol or chlorine, as this could damage aluminium and polyamide materials.
- Dry with a clean, dry soft cloth.
- Make sure there is no accumulation of dirt or debris around the axle, brakes or wheels.
- Periodically check all screws and bolts for wear and strength. Replace all broken, damaged or worn parts immediately.
- In the event of a malfunction, do not repair the table yourself. Consult an authorised dealer.

CHANGE OF OWNERSHIP

The product can be reused. It must first be cleaned and disinfected in accordance with hygiene instructions. When transferring the rollator, remember to hand over all the necessary technical documents to the new user. Beforehand, the product must be inspected by an authorised specialist.

RECYCLING

When the product has become unusable and you have to throw it away: please contact your specialist dealer. If you wish to carry out recycling for yourself, ask your local waste management company for disposal guidelines.

STORAGE

Failure to comply the storage conditions may lead to deterioration of the product and thus to the risk of serious injury.

Do not store your product for an extended period of time near a heat source or in the sun (e.g. behind a window or near a radiator) or near a cold source.

Keep away from all flames and sources of sparks.

Respect the conditions for storing and storing your product.

- In a dry and temperate place
- Protect your product by packaging it from dust, corrosion (e.g. abrasive elements, sand, sea water, salty air).
- Store all removed parts together in the same place (or locate them if necessary) to avoid mixing them with other products when reassembling.
- All components must be stored without any load (do not place too heavy objects on the product components, do not jam anything...).

TECHNICAL SPECIFICATIONS

| | |
|----------------------------------------|--------------------------------------------|
| Overall dimensions | Length 82 x width 45 x height 79 to 114 cm |
| Height ajustable | 78 cm to 105 cm |
| Dimensions of the large tray | 39 cm x 80 cm |
| Passage between the U-shaped structure | 59 cm |
| Total weight | 13 kg |
| Maximum supported weight | 25 kg |

PERIOD OF USE

The life expectancy of this product is 5 years under normal use, safety and maintenance conditions. Beyond this period, the product can be used for as long as it is still in good condition.

WARRANTY

This product is guaranteed for two years from the date of purchase by the user. IDENTITIES will provide replacement parts for any part that has broken due to a manufacturing, design or material defect. Contact the distributor from whom you purchased the product in case you need to make a claim. Failures resulting from accidental damage, misuse, alteration, or wear and tear are not covered by this warranty. Contact the dealer from whom you purchased the product to apply the warranty.

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